

# Rabbit Monoclonal Primary Antibodies

Helping pathologists with tools for a timely and accurate diagnosis

HER-2/neu Rb Mono

We now provide


**PATHWAY<sup>®</sup> HER-2/neu**

**HER-2/neu (4B5)**






**FDA Approved for VIAS<sup>™</sup>**

(Ventana Image Analysis System)

FDA Approved

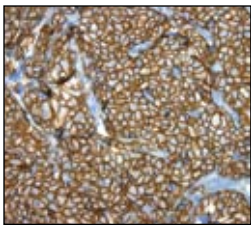

HER-2/neu (4B5)

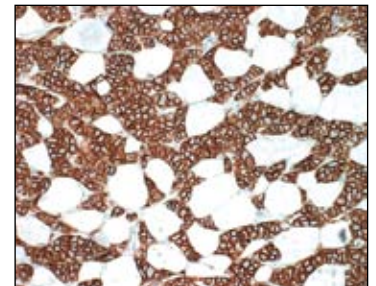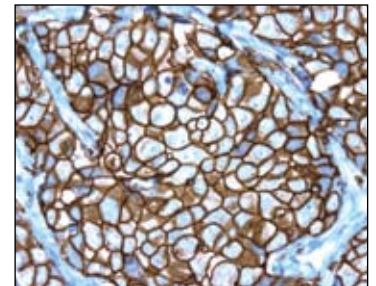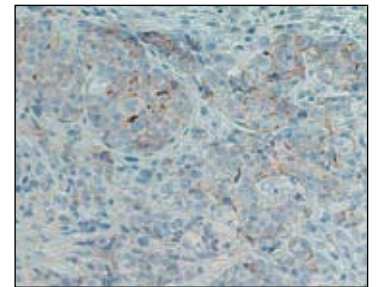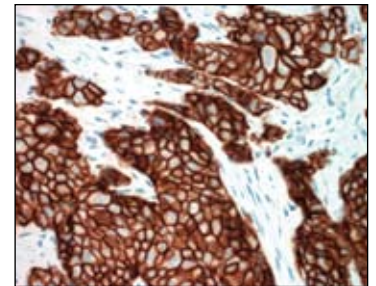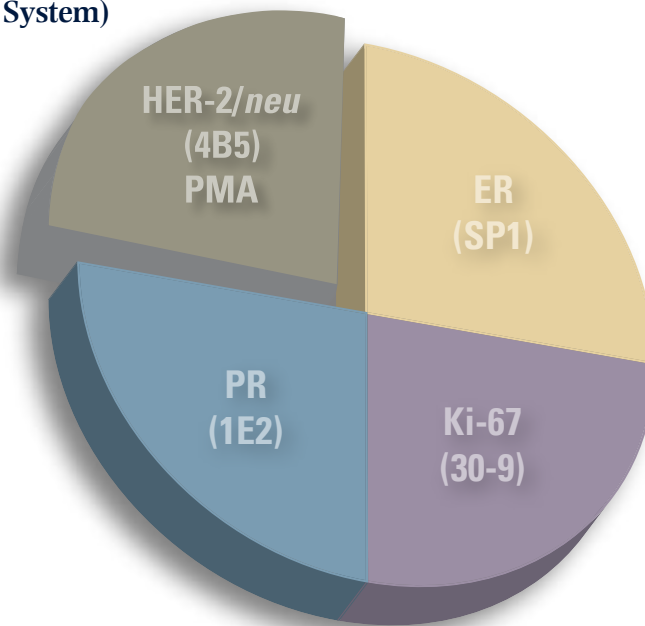

## The Advanced Complete Choice


in Moving from Mouse to Rabbit Monoclonal Primary Antibodies for Breast Cancer

- Ventana Medical Systems' (Ventana) PATHWAY anti-HER-2/neu (4B5) Rabbit Monoclonal Primary Antibody is a rabbit monoclonal antibody intended for laboratory use for the semi-quantitative detection of HER2 antigen in sections of formalin-fixed, paraffin-embedded normal and neoplastic tissue on a Ventana automated immunohistochemistry slide staining device.
- It is indicated as an aid in the assessment of breast cancer patients for whom Herceptin<sup>®</sup> treatment is considered.
- PATHWAY anti-Her-2/neu is a rabbit monoclonal primary antibody (clone 4B5) directed against the internal domain of the c-erbB-2 oncoprotein (HER2).
- This antibody has been optimally diluted for use with Ventana iView<sup>™</sup> DAB and ultraView<sup>™</sup> detection kits.
- This antibody is intended for *in vitro* diagnostic (IVD) use.



## PATHWAY HER-2/*neu* (4B5) Rabbit Monoclonal Primary Antibody

Type of Test	Immunohistochemistry Antibody
Part Numbers	790-100 Kit (U.S. ONLY) <i>Includes: PATHWAY HER-2 4in1 control slides (10), Interpretation Guide 790-2991 (50 tests) Outside U.S.</i>
Clone	Rabbit Monoclonal 4B5
Quantity	5 mL dispenser (50 tests)
Storage	2°- 8° C
Run Time	16-32 minutes primary antibody incubation on NexES <sup>®</sup> and BenchMark <sup>®</sup> Series
Regulatory Status	IVD PMA

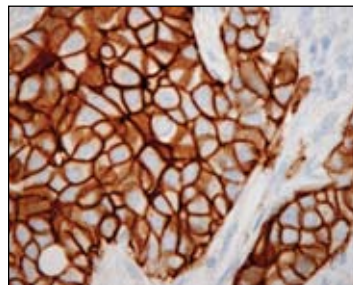

### Experience Other Ventana Products:



<i>ultra</i> VIEW <sup>™</sup> Universal DAB Detection Kit	Catalog No. 760-500
<i>i</i> View <sup>™</sup> DAB Detection Kit	Catalog No. 760-091
Hematoxylin II	Catalog No. 790-2208

## The Complete Choice

Ventana Medical Systems offers a comprehensive selection of quality IHC, ISH, and FISH reagents, as well as automated slide processing capabilities.

Automation for Consistent and Reproducible Results

Extensive Selection of Antibodies, Probes and Detection Reagents

Flexible System Software for Creating Custom Assay Protocols



### Corporate Offices

**North America**  
Ventana Medical Systems, Inc.  
1910 E. Innovation Park Drive  
Tucson, Arizona 85755  
U.S.A.  
+1 (520) 887-2155  
(800) 227-2155 (U.S.)

**Europe**  
Ventana Medical Systems, S.A.  
Parc d'Innovation - BP 30144  
Rue G. de Kaysersberg  
F - 67404 Illkirch CEDEX  
France  
+33 (0) 3 90 40 52 00

**Japan**  
Ventana Japan K.K.  
ベントナ・ジャパン株式会社  
〒220-8135 神奈川県横浜市  
西区みなとみらい 2-2-1  
横浜ランドマークタワー35 階  
+81 (0) 45 228 5071

**Australia, New Zealand**  
Ventana Medical Systems Pty Ltd  
5/39 Grand Boulevard  
Montmorency VIC 3094  
Australia  
+61 (0) 3 9431 6064

EC	REP
----	-----

**MDCI Ltd.**  
Arundel House  
1 Liverpool Gardens  
Worthing  
West Sussex BN11 1SL  
U.K.